

Reduce waste. Go green.


The Cardinal Health reprocessing solution is designed to meet sustainability goals and reduce waste.

Kendall SCD™ Sequential Compression Sleeves are intended for single-patient use. Used compression sleeves should be placed in the provided collection hampers after the patient is discharged.

After patient use:



Detach tubing from the Kendall SCD™ Sequential Compression Sleeves.

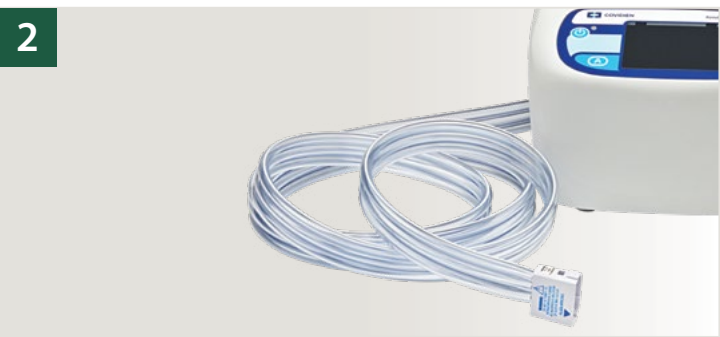

Leave tubing connected to the Kendall SCD™ controller.



Deposit the Kendall SCD™ Sequential Compression Sleeves in the hamper.

For more information call **844.268.2583**
or email **gogreen@cardinalhealth.com**



Do not use collection hampers to discard the following items



Sharps
No scalpels, needles, syringes, sutures or glass



Waste
No personal equipment, trash, patient care disposables or linens



Batteries
Remove from device but do not place in collection container



Tubing or Kendall SCD™ controllers
Leave tubing connected to the Kendall SCD™ controllers



**Sustainable
Technologies™**
A Cardinal Health business