

Mail-back Collection Bin Service Guide

This document will review the processes involved in setting up, packaging and shipping various collection containers utilized in the Sustainable Technologies™ Mail-back Program.



Hamper with green collection bag

Container setup



Hands-free lid

Pedal

Used Medical Device
Caution: For processing, handling or incineration only

Sustainable Technologies™
A Cardinal Health business

Sustainable Technologies, 8322 Eagle Point Drive, Suite 410, Riverside, FL 33516
For service please call: (844) 208-2563

- 1 Press pedal
- 2 Fold bag over inner frame
- 3 Pull edges down four inches



Container packaging and shipping



Sustainable Technologies™
A Cardinal Health business

Sustainable Technologies, 8322 Eagle Point Drive, Suite 410, Riverside, FL 33516
For service please call: (844) 208-2563

- 1 Tie bag with twist tie
- 2 Place bag into corrugate shipping box supplied by Sustainable Technologies
- 3 Ship back to Sustainable Technologies using supplied shipping label



2 gallon wall-mounted container, 2 or 5 gallon container

Container setup

2 gallon wall-mounted container



1 Mount bracket on wall

2 Place container in bracket



2 or 5 gallon container



1 Fix back of lid to container

2 Place container in collection location

Container packaging and shipping



1 Place provided lid with shipping label on collection container




2 Seal container lid and ship

Blue 18 gallon container with wire cart

Container setup



Container packaging and shipping

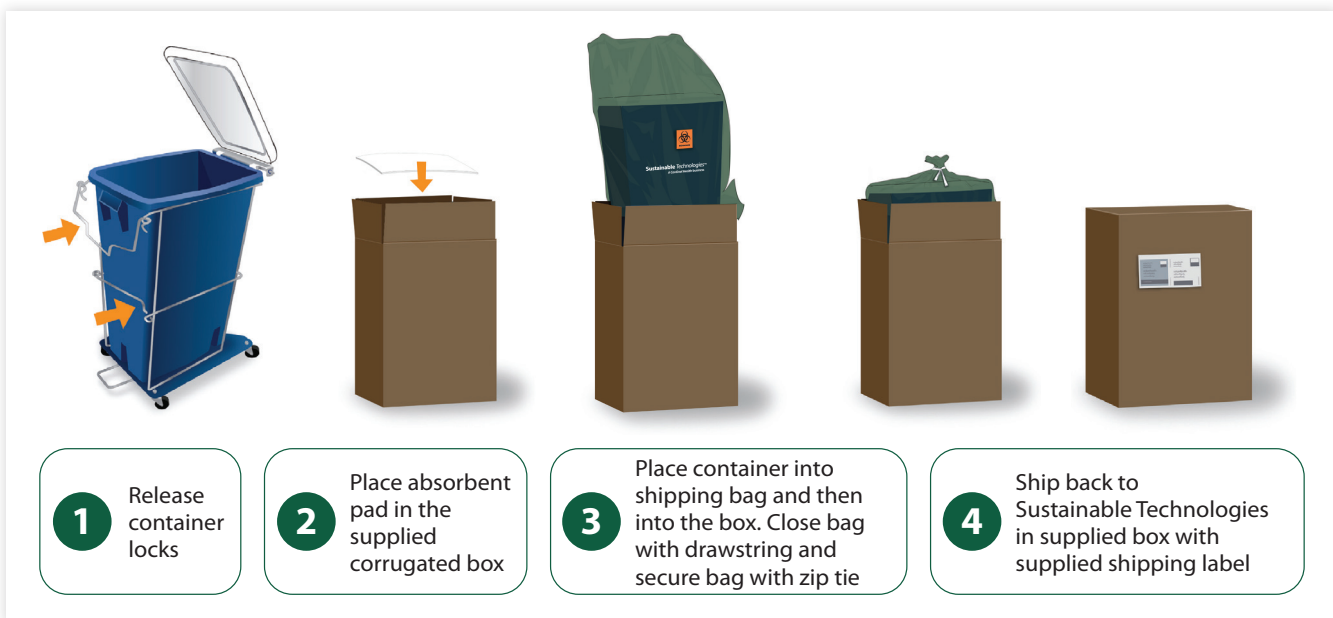

For more information, call **844.268.2583** or email **gogreen@cardinalhealth.com**

Sustainable Technologies™
A Cardinal Health business

sustainabletechnologies.com

© 2023 Cardinal Health. All Rights Reserved. CARDINAL HEALTH, SUSTAINABLE TECHNOLOGIES and the Sustainable Technologies LOGO are trademarks of Cardinal Health and may be registered in the US and/or in other countries. All other marks are the property of their respective owners. Lit no. 25T23-2701167 (12/2023)

7000 Cardinal Place
Dublin, Ohio 43017